



Guardian

The electronic witnessing system using dual identification label technologies, designed for seamless protection for patients and IVF clinics.



[Patient list](#)
[Batch Printing](#)
[Comprehensive Inventory](#)
[Statistical Analysis](#)
[Consumables Management](#)
[System Configuration](#)

Please select ↕ ↻

[Full Process Mode](#)
[Console Mode](#)

Active Cycle

Remaining 457 cycles

Incoming	Day 0	Day 1	Day 2	Day 3	Day 4	Day 5	Day 6	Day 7
Cook Virginia 1217 Fresh Cycle, 1.11	Kingsley Bronwyn 1207 Fresh Cycle, 11.20	Ball Lucille 1185 Thaw Cycle, 3.12	Duck Daisy 1042 Fresh Cycle, 1.4	Ball Lucille 1185 Thaw Cycle, 3.15	Brown Corinne 1099 Thaw Cycle, 11.5	Carr Alicia 12300 Thaw Cycle, 2.3	Doe Jane 9285 Thaw Cycle, 4.1	Edwards Katrina 0962 Thaw Cycle, 10.1
McDonald Sydney 4164 Fresh Cycle, 1.12	Blackburn Jasna 1264 Fresh Cycle, 11.23	Anderson Apple 0194 Thaw Cycle, 3.12	Blackburn Jasna 1264 Fresh Cycle, 11.3	Albright Sally 1170 Thaw Cycle, 3.3	Cisak Amelia 1142 Thaw Cycle, 2.2	Estevéz Maria 9866 Thaw Cycle, 10.2		
Duck Daisy 1042 Fresh Cycle, 1.10	Albright Sally 1170 Thaw Cycle, 3.14	Cook Lily 61202 Fresh Cycle, 2024090.4	Anderson Apple 0194 Thaw Cycle, 3.2		Casis Belicia 1189 Thaw Cycle, 11.10			
Olsen Sandy 1028 Fresh Cycle, 11.12	Betting Becky 0857 Thaw Cycle, 3.11				Brown Melissa 2546 Thaw Cycle, 2.1			

INTRODUCING GUARDIAN

Guardian is an electronic witnessing system using both types of identification label technologies, RFID chip and QR codes, designed as a seamless protection system for patients and IVF clinics. The customer can choose the label options for each ID check procedure to achieve an optimised workflow. The hand-held PDA can read both types of labels which allows for ID checks to be done anywhere and everywhere. Guardian aims to help labs manage daily workflow and ensure ID checks are performed per protocol. Guardian is your reliable lab safety management assistant.



Customisable System

The one and only system on the market that allows customers to use RFID or QR codes flexibly.



Seamless Integration


Hardware has a small footprint, is portable, and enables seamless implementation into labs.



Diverse Workflow Display

ID-check status and workflow can be viewed from D0 to D7.

GUARDIAN: THE INNOVATIVE TRACKING TECHNOLOGY LEADER

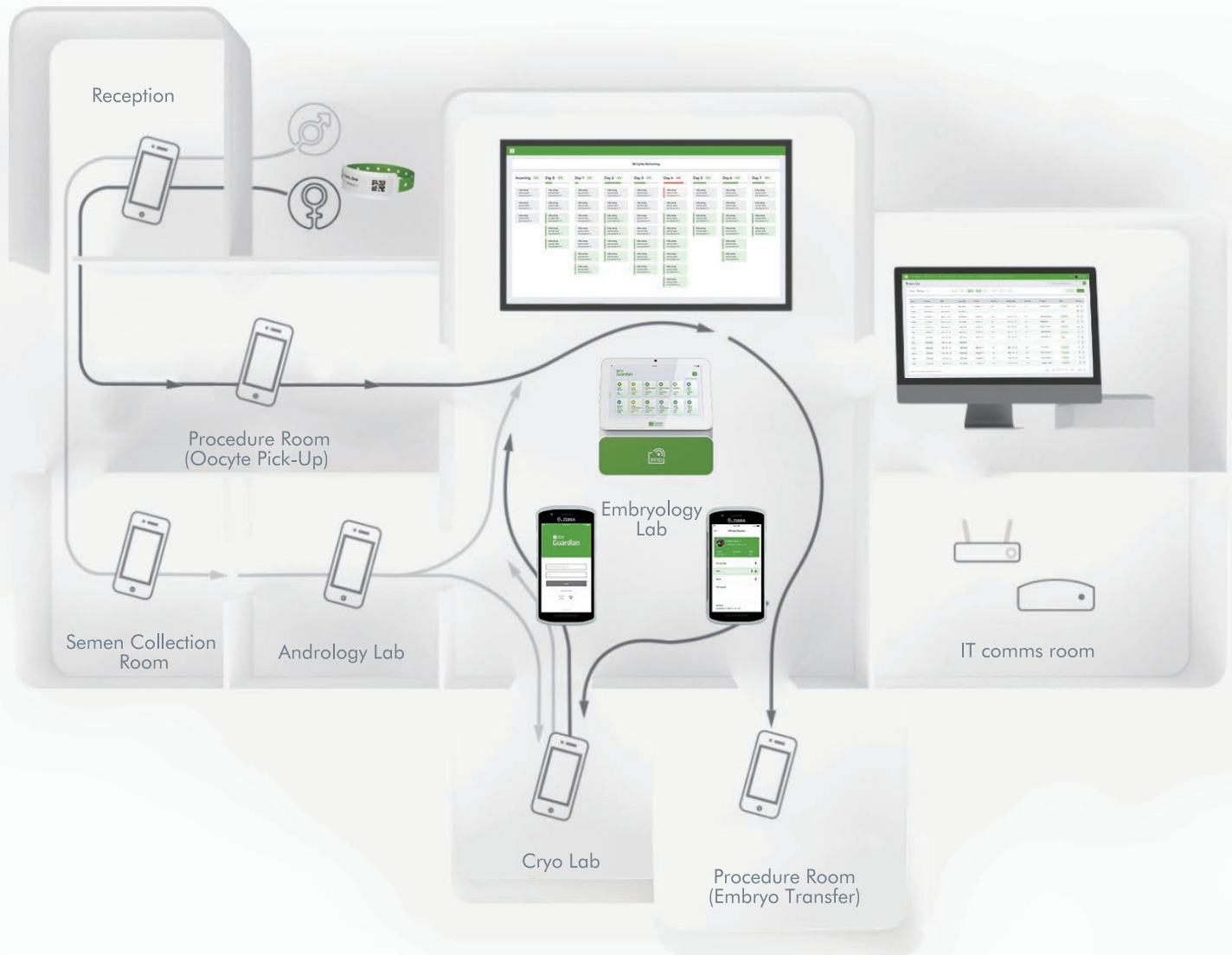



Personal Digital Assistant (PDA)

Witnessing on the move for RFID and QR codes.



Workstation

Guardian can be used in desktop mode in a network connected environment.



RFID All-in-one

Touch screen workstation with RFID reading inbuilt.



Monitor

Project guardian to an external lab display to quickly monitor progress and alarms.



Label Printer

Enabling printing of RFID tags, QR, and cryogenic labels for your workflow.



Server

The beating heart of guardian containing software & databases, with connectivity to the clinic network.

GUARDIAN: VERIFIES PATIENT ID ACROSS ALL STEPS AND LOCATIONS



Advanced Safety


Write electronic information and print visual labels simultaneously on RFID tags.

Flexibility Delivered

Select RFID or QR codes where needed.

Undisruptive Installation

No physical hardware change to existing equipment and lab set up.

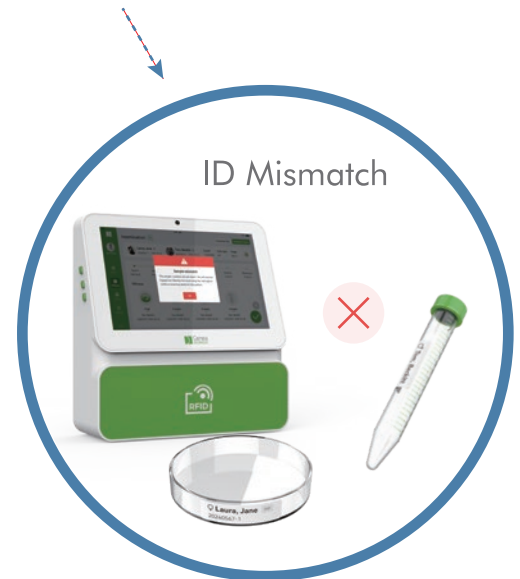


Mobility Built-In

Mobile for witnessing to be done anywhere and everywhere.



Transparency Connected

Mismatches distributed in All in one, PDA and dashboard for full visibility.

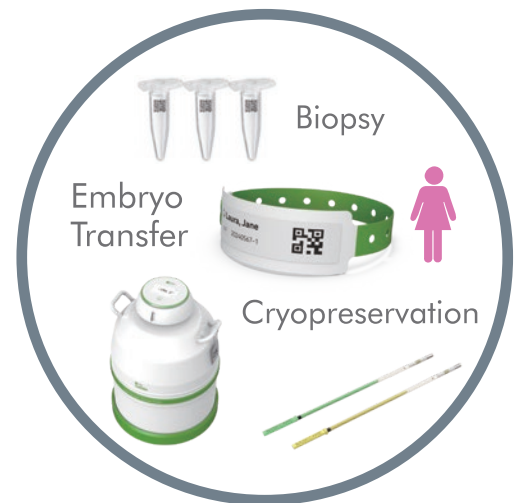

For healthcare professionals only. Please refer to the instructions for use.
For further information, please contact your Service Representative or visit:
www.geneabiomedx.com