



Guardian

The electronic witnessing system using dual identification label technologies, designed for simplified sample traceability for patients and IVF clinics.



Admin Set | Batch Printing | Comprehensive Directory | Statistical Analysis | Consumables Management | System Configuration

Please select

Jason Taylor Embryologist

Active Cycle

Full Process Mode Console Mode

Remaining 457 cycles

Incoming	Day 0	Day 1	Day 2	Day 3	Day 4	Day 5	Day 6	Day 7
Cash Virginia 1202 Thaw Cycle, 1, 1	Kingsley Bronwyn 1187 Thaw Cycle, 11, 20	Bull Lucille 1185 Thaw Cycle, 1, 11	Duck Daisy 1041 Thaw Cycle, 1, 4	Bull Lucille 1185 Thaw Cycle, 1, 5	Brown Corinne 1099 Thaw Cycle, 11, 3	Carr Alicia 12300 Thaw Cycle, 2, 1	Doe Jane 9283 Thaw Cycle, 4, 1	Edwards Katrina 9902 Thaw Cycle, 20, 1
McDonald Sydney 492 Thaw Cycle, 1, 2	Blackburn Jaana 8194 Thaw Cycle, 11, 21	Anderson Apple 8194 Thaw Cycle, 3, 12	Blackburn Jaana 8194 Thaw Cycle, 11, 8	ABright Sally 2176 Thaw Cycle, 3, 1	Cole Amelia 1242 Thaw Cycle, 2, 2			Estevez Maria 9686 Thaw Cycle, 10, 2
Duck Daisy 1041 Thaw Cycle, 1, 10		ABright Sally 2176 Thaw Cycle, 3, 14	Cook Lily 14081 Thaw Cycle, 20, 10, 1	Anderson Apple 4124 Thaw Cycle, 3, 2		Casin Belcia 1188 Thaw Cycle, 11, 10		
Olsen Sandy 828 Thaw Cycle, 11, 2		Betting Becky 887 Thaw Cycle, 2, 11				Brown Melissa 2148 Thaw Cycle, 2, 1		



INTRODUCING GUARDIAN

Guardian is an electronic witnessing system using both types of identification label technologies, RFID chip and QR codes, designed as a simplified sample traceability system for patients and IVF clinics. The customer can choose the label options for each ID check procedure to achieve an optimised workflow. The hand-held PDA can read both types of labels which allows for ID checks to be done anywhere and everywhere. Guardian aims to help labs manage daily workflow and ensure ID checks are performed per protocol. Guardian, offering seamless traceability everywhere you need it.



Customisable System

It seamlessly adapts to the workflow of your lab, rather than requiring your lab to adapt to it.



Seamless Integration


Hardware has a small footprint, is portable, and enables seamless implementation.



Diverse Workflow Display

ID-check status and workflow can be viewed from day zero to day seven.

GUARDIAN: SEAMLESS TRACEABILITY, EVERYWHERE YOU NEED IT

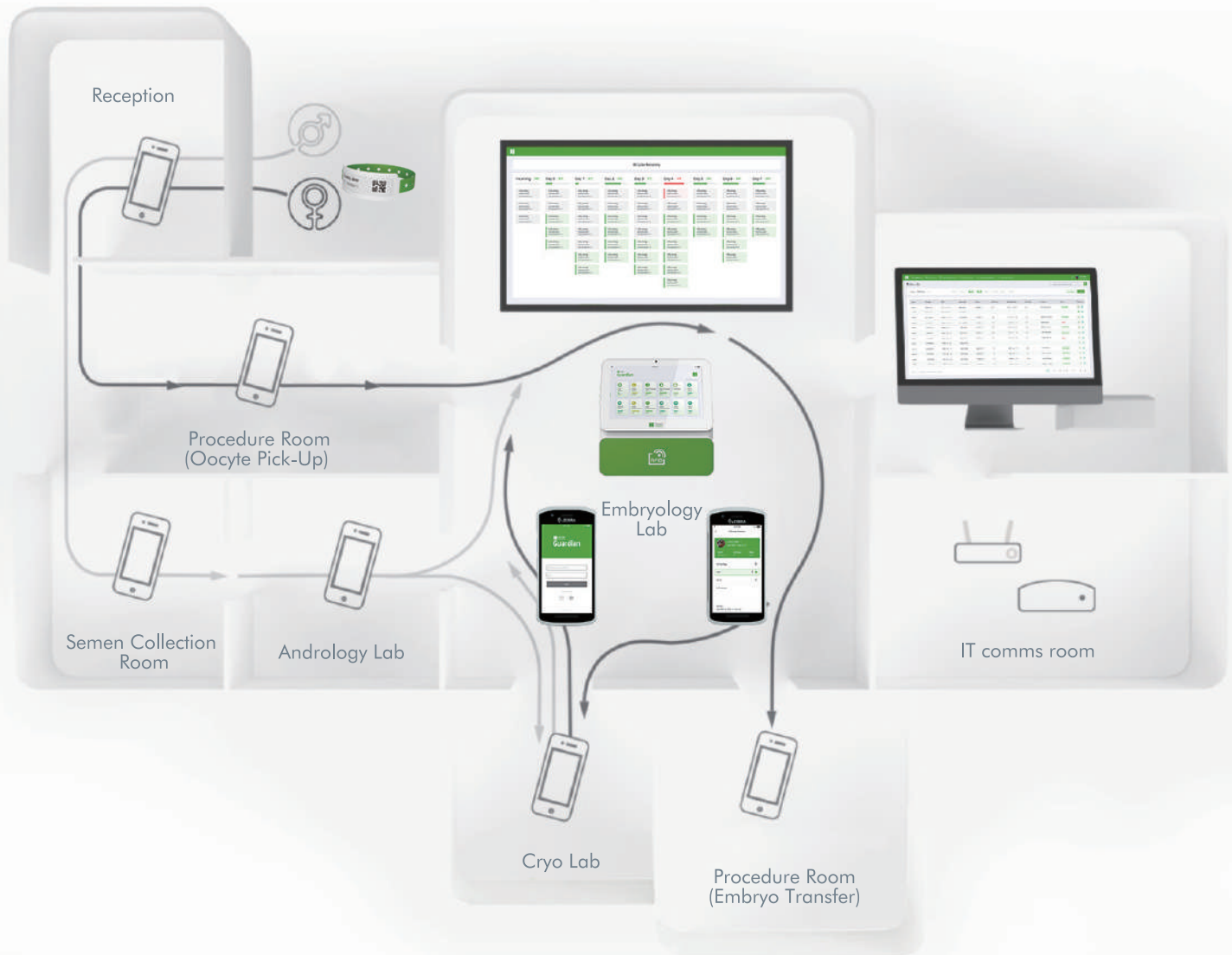



Personal Digital Assistant (PDA)

Witnessing on the move for RFID and QR codes.



Workstation

Guardian can be used in desktop mode in a network connected environment.



RFID All-in-one

Touch screen workstation with RFID reading built-in.



Monitor

Project guardian to an external lab display to quickly monitor progress and alarms.



Label Printer

Enabling printing of RFID tags, QR, and cryogenic labels for your workflow.



Server

The beating heart of guardian containing software & databases, with connectivity to the clinic network.

GUARDIAN: VERIFIES PATIENT ID ACROSS ALL STEPS AND LOCATIONS



Seamless Traceability


Write electronic information and print visual labels simultaneously on RFID tags.

Flexibility Delivered

Select RFID or QR codes where needed.

Undisruptive Installation

No physical hardware change to existing equipment and lab set up.

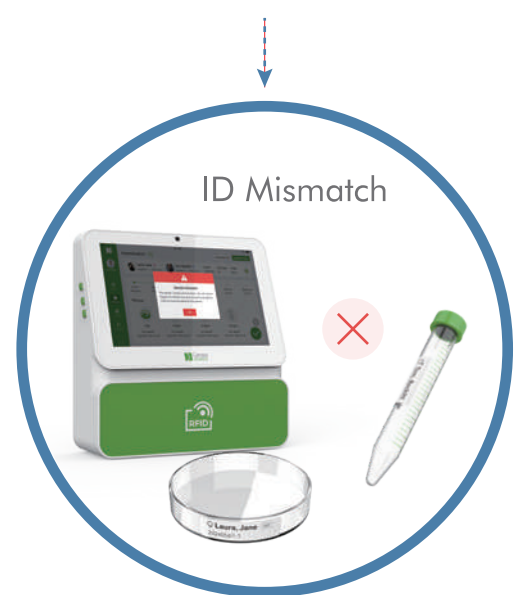


Mobility Built-in

Mobile witnessing anytime and everywhere.



Transparency Connected

Mismatches distributed to PDA and dashboard for full visibility.

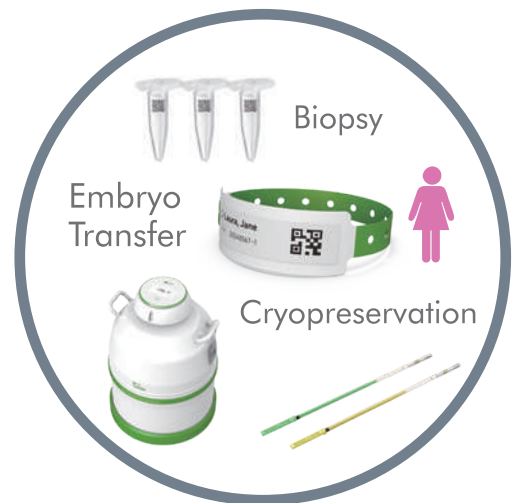

For healthcare professionals only. Please refer to the instructions for use.
For further information, please contact your Service Representative or visit:
www.geneabiomedx.com